

EASYFLEX N NON-PVC IV BAG

Reduces the risk of needle stick injuries

Allows for the preparation and administration of intravenous therapies in an easy and safe way

With a needle-free injection port and a universal infusion port




spiritmedical
supplying for life

GET IN TOUCH TO FIND OUT MORE ABOUT OUR RANGE OF IV FLUID BAGS

✓ **SAFE**

Needle-free injection ports reduce the risk of needle stick injuries

✓ **QUICK**

Allows for faster transfer of solutions & medicines

✓ **EASY**

Easy to use and compatible with luer lock connectors and syringes

EasyFlex N Non-PVC IV Bag

A clear, flexible needle-free IV fluid bag from our safe infusion range, which allows for the preparation and administration of intravenous therapies in an easy and safe way. It features an integral needle-free access port and a universal infusion port.

KEY FEATURES & BENEFITS

- ⊕ Clear see through bag
- ⊕ Flexible and lightweight
- ⊕ PVC free (made from polyolefin)
- ⊕ Latex free
- ⊕ DEHP free
- ⊕ Leak free
- ⊕ Single use sterile bag
- ⊕ Simple, quick and easy to use
- ⊕ One needle-free injection port
 - Compatible with all luer lock connectors
 - Compatible with luer lock syringes
 - Simply push and twist for secure connection
- ⊕ One universally compatible infusion port
 - Simply twist off the protective cover and fully insert the giving set




Hanging eyelet

Film: Polyolefin
(Ethylene/Propylene Copolymer)

Infusion port
twist-off cap

Injection Port
Needle-free access
Polycarbonate
and silicone

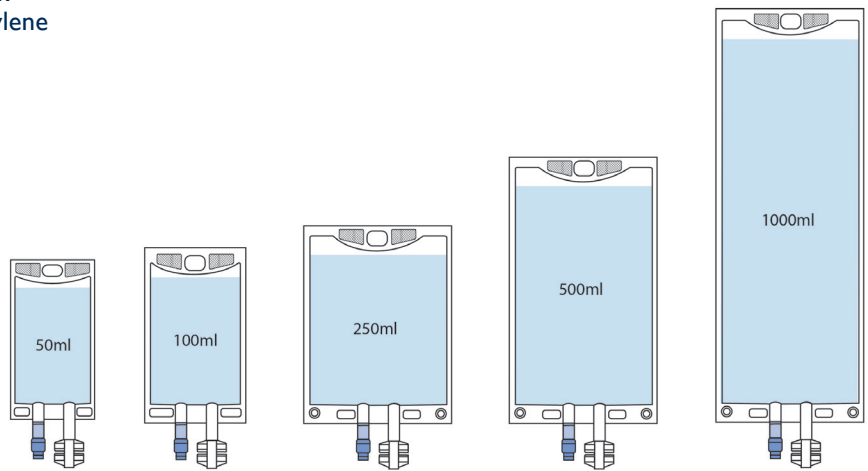

SOLUTION	VOLUME	PACK QUANTITY	PRODUCT CODE
Sodium Chloride 0.9%	50ml	70 bags	ACG0116GB
Sodium Chloride 0.9%	100ml	60 bags	ACG0216GB
Sodium Chloride 0.9%	250ml	30 bags	ACG0416GB
Sodium Chloride 0.9%	500ml	18 bags	ACG0716GB
Sodium Chloride 0.9%	1000ml	10 bags	ACG1016GB
Glucose 5%	50ml	70 bags	AFG0116GB
Glucose 5%	100ml	60 bags	AFG0216GB
Glucose 5%	250ml	30 bags	AFG0416GB
Glucose 5%	500ml	18 bags	AFG0716GB

TO ORDER OR FOR FURTHER INFORMATION PLEASE VISIT OUR WEBSITE